

## SKINTACT Cool Contact • Electrosurgical Universal Neutral Electrodes


“SKINTACT Cool Contact Electrosurgical Neutral Electrodes provide safety of the highest degree. They are compatible with all common generators on the market and fully comply with relevant clauses of EN 60601-2-2:2009 guaranteeing patient safety and your peace of mind.”



# SKINTACT Cool Contact • Electrosurgical Universal Neutral Electrodes

## The professional's choice

- ❑ round leading edge geometry disperses current more efficiently thus maintaining a cooler interface between patient and neutral electrode
- ❑ overlapping gel reduces risk of incident as aluminum edge exposure is eliminated
- ❑ softer hydrogel is kinder to skin and creates a better interface between patient and neutral electrode facilitating a more efficient electrical contact

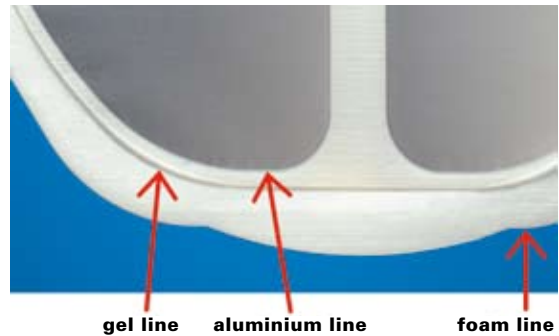

## SKINTACT Cool Contact Electrosurgical Universal Neutral Electrodes



Order Code	RO01	RO21
Type	unsplit	split
Size	universal	universal
Conductive Area	118 cm <sup>2</sup>	107 cm <sup>2</sup>
Packaging Standard	1 or 5 / pouch 50 / box	1 or 5 / pouch 50 / box

with 3 m pre-attached disposable cable*		
Order Code	RO01B30	RO21A30
Connector	Non-REM	REM compatible
Packaging Standard	1 / pouch 50 / box	1 / pouch 50 / box

**Compatibility:** In any case, split Skintact Neutral Electrodes are compatible with the dynamic contact quality monitoring systems REM™, NESSY™ and ARM™.

- ❑ universal for use on pediatric and adult
- ❑ user friendly tracking system (removable labels)
- ❑ shoulder tab to enable theatre staff to easily lift the neutral electrode off its protective liner without damaging gel or creasing aluminum
- ❑ bottom end tab to facilitate removal from patient after procedure with minimal skin trauma



**SKINTACT Cool Contact Reusable Dispersive Cable with stamped Lot Number for easy tracking (types REM and NON-REM).**

## SKINTACT Cool Contact Reusable Dispersive Cables\* and Adaptors

Order Code	NEK300A	NEK300B	NEK300D	NEK300E	NEK300G	NEK300H
Connector (compatible)	REM (Valleylab)	Non-REM (Valleylab)	6.3 mm jack (standard)	4 mm banana	2 Pin adjustable (Eschmann)	6.3 mm jack (Erbe)
Cable	3 m	3 m	3 m	3 m	3 m	3 m

\* SKINTACT Cables can also be supplied in 5 metre.

Order Code		NEKA-AD		NEKA-AH		RAZOR
Description		Adaptor A/B-D		Adaptor A-H		Prep razor

SKINTACT Cool Contact Electrosurgical Neutral Electrodes fully comply with relevant clauses of EN 60601-2-2:2009.